

## **HairMax LaserComb**

### **A Typical User Experience.**

As we all know, skin cells need light. Your hair is no different. The powerful, safe and stimulating laser light in HairMax LaserComb gives energy to your scalp and hair follicles.

Through the scientific principles of PhotoBioStimulation, Laser Phototherapy progressively improves the quality, strength and thickness of your hair over time. Laser light therapy, also called photo therapy in the medical community, is a safe, effective therapy with no known side effects.

Indeed, international experiences and studies of the efficacies of Laser Phototherapy **go far beyond** those that are described on our web site.

Results vary from individual to individual, but within the first 6 to 12 weeks, you should notice early general improvements or activation of your hair occurring. These improvements include a normalizing of scalp condition, shinier and thicker hair. Then, over the next few months, you may notice your hair gradually becoming fuller and stronger. Over time, all users – men and women – usually notice progressively increasing benefits as they use their HairMax LaserComb. Conversely, once the usage stops, our experience is that hair and scalp conditions gradually regress back to where they were.

Like brushing your teeth to maintain healthy teeth & gums, regular HairMax usage is recommended. As little as once or twice a week for 10 minutes may be sufficient to maintain achieved results over time.

## IMPORTANT INSTRUCTIONS ON USING HAIRMAX


**We recommend using the LaserComb for a minimum of 5 minutes or a maximum of 15 minutes per treatment**

- 3 times a week for the first 4 weeks
- 2 times a week thereafter

It is recommended that the LaserComb be used on a clean scalp after washing. Your hair can be wet, towel dried or dry. The LaserComb should be used before applying any other topical treatments, gels, sprays or concealing products..




1. Insert the power attachment into the LaserComb.

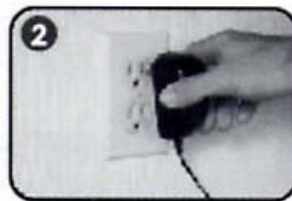

2. Plug in the wall adaptor.



3. Turn on the LaserComb by first switching on the 'Power' button followed by the 'Laser' button.

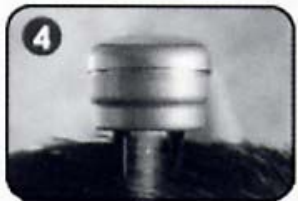

4. Place the LaserComb flat on your scalp using the two rows of teeth.



5. Move the LaserComb using a 'Step' method by placing the LaserComb on a spot and leaving it there for 4 or 5 seconds then moving it 1/2 inch to the next spot. (Mentally counting "one thousand, two thousand, etc." should help with timing).

- 6 Move the LaserComb against the direction of the hair growth. With long hair, use an ordinary comb or your hand to lift the hair in advance of the LaserComb. Make 2 complete passes with the LaserComb over the whole scalp area in the manner described above.
- 7 After treatment turn off the power button, remove the wall plug and return the LaserComb and power supply to the case for storage.

### Please Note

- The total time treatment with the LaserComb should take will range between 5 to 15 minutes.

- It is important to move the LaserComb very slowly. If you move it too fast, the treatment will not be as effective.
- Like any therapy, it is recommended that you do not overuse the LaserComb as additional usage may be less effective.


**IMPORTANT SAFETY INSTRUCTIONS**  
**-READ CAREFULLY BEFORE USE-**

**LASER SAFETY** Warning and Cautions:

1. Do not stare directly at the laser light or its reflection off a mirror-like surface (as this can act like a direct beam on the eye).
2. Do not shine the laser light in or around your eyes
3. Keep out of reach of children. Do not allow minors to use lasers.
4. Never shine a laser at anyone.
5. The use of magnifying optical instruments with this product can cause eye danger if the laser is shined directly into the eye.
6. Do not attempt to repair, enter, alter or interfere in any way with the workings of this appliance or any part of this appliance.
7. Only use directions, cleaning methods and operating procedures specified in this manual.

**Note:** This appliance has no user-serviceable parts.

**The following laser safety label can be found on this product:**

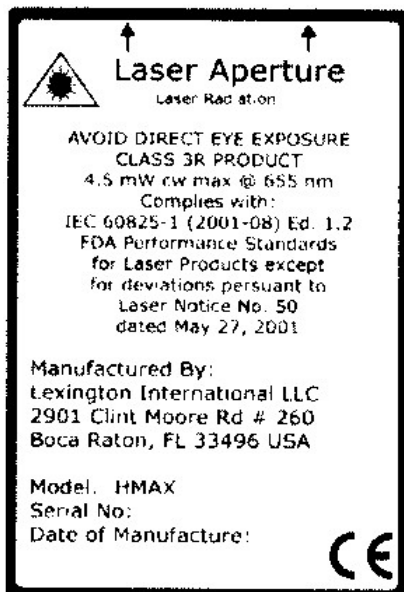

Location: On the underside of the handle

## **ELECTRICAL SAFETY**

Basic safety precautions should be observed when using this, or any other electrical device, especially when children are present.

**DANGER:** Keep away from water. As with most electrical appliances, contact with water may result in an electrical shock.

To reduce the risk of hazard by electric shock.

- Always unplug the appliance from the power outlet immediately after using.
- Do not use while bathing.
- Do not place or store this appliance where it can fall or be pulled into a tub or sink.
- Do not expose this appliance to running water, drop into water or any other liquid.
- If the appliance falls into water, unplug it immediately from the power outlet. Do NOT reach into or touch the water.
- If using this appliance near water, the use of an electrical outlet protected by a Ground Fault Circuit interrupter is recommended.

**WARNING:** To reduce the risk of burns, electrocution, fire or injury:

- This appliance should never be left unattended when plugged into a power outlet.
- Care should be taken when operating this appliance near children. Children are naturally curious and this appliance is not a toy and should not be played with by minors.
- Use this appliance only for its intended use as described in this manual.
- Never operate this appliance if it has a damaged cord or plug. If it is not working properly, or if it has been dropped, damaged or immersed in water, return the appliance to Lexington International, LLC for examination, service or replacement.
- Keep the cord away from heated surfaces.

**NOTICE:** If any malfunction occurs, do not attempt to repair it yourself.

This appliance has no user-serviceable parts.

### **Power Supply Specifications:**

Output from transformer: 9V, minimum 300mA. Plug: 2.1mm x 5.5mm.

For use outside TW, please make sure the power supply adapter matches these specifications.

## MAINTENANCE INSTRUCTIONS


- Always unplug this appliance from the power outlet while cleaning it.
- Warning! Do NOT use acetone or any other solvents on any part of this appliance. Use of acetone or other solvents to clean the window can immediately block most of the laser light from coming through.
- To clean the laser window (from where the laser light emanates), gently clean with a soft, slightly moistened cloth.
- If the cord becomes twisted, unplug and straighten it prior to use.
- Do not put strain on cord if you need to untangle or straighten it.

### **Cleaning comb teeth:**

To keep the teeth of the combs clean, gently remove any hairs or dust with a slightly moistened towel.

### **Multi-user hygiene:**

For use by more than one person, a second set of combs have been included in the carrying case. Gently pull the combs down to install the second set of combs.

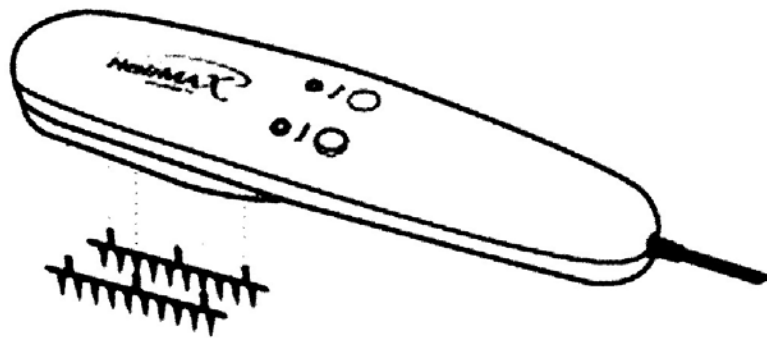

**For a more thorough cleaning, gently pull out hair combs and replace them using a little gentle manipulation.**

---

**SAVE THESE INSTRUCTIONS**

---